




For All Your Clinical Diagnostic Needs

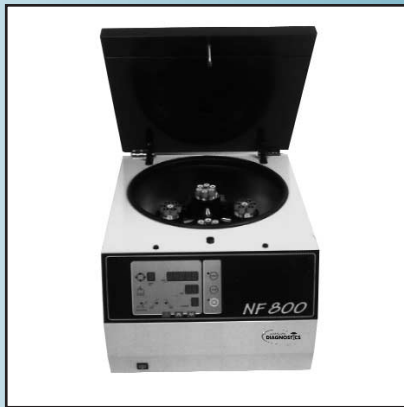

*Capable of handling all Benchtop applications • Durable
Transparent lid, safety interlocks, no exposed parts, sealed dome shields
Runs only a few degrees higher than ambient temperature • Can hold six tubes*

**Digital fully programmable, Haematocrit, and Multi-Purpose
benchtop centrifuges also available**

BENCHTOP CENTRIFUGE

SYSTEM FEATURES:

Capable of handling all Benchtop applications. It is virtually indestructible. The Benchfuge was designed for:

- Physicians offices
- Teaching labs
- Small Laboratories
- R & D laboratories
- Insurance examiners
- Satellite labs

Built of high impact ABS and Lexan, to survive the rigors of daily use. Lexan lid with interlock and sealed dome shields minimize biohazard risks.

No exposed moving parts. Lexan tube holders are easily removed and sterilized. A self-contained one piece head prevents tubes from falling through.

Accommodates 6 standard 10ml to 15ml (16 x 75 mm and 16 x 125 mm) tubes as well as smaller sizes with optional tube inserts

Inserts can be sterilized in case of breakage of a tube

Optional electronic timer model to set exact spin times. No guess work as with mechanical timers

Optional fully automatic timer brake system provides smooth full stopping in under 30 seconds

Benchfuge assures operating temperatures only a few degrees higher than ambient temperature no matter how long the centrifuge runs.

TECHNICAL SPECIFICATIONS:

- The Benchfuge conforms to all UL, FDA, OSHA, and VDA regulations
- Power requirements 115 VAC, 60 HZ, 1.5 amps or 220 VAC, 50 HZ, .75 amps
- Rated Speed (RPM- Max): 3400
- Relative Centrifuge Force: 1450 nominal
- Cord: UL listed, 6-ft. (1.8m), strain-relieved, grounded cord and plug
- Dimensions: 11 x 11 inches (28 x 28 cm)
9 _ inches high (23 _ cm)
- Number of places: 6
- Angle of rotation: 49(deg)
- Weight: 7lbs. (3.2 kg)
- Shipping weight: 10 lbs. (4.5 kg)
- One year warranty

PLEASE CONTACT FOR SPECIFICATIONS OF HAEMATOCRIT, DIGITAL FULLY PROGRAMABLE, AND MULTI-PURPOSE CENTRIFUGES

PLEASE CONTACT:

WORLDWIDE DIAGNOSTICS

P.O. Box 656 • Ridgefield, CT 06877 • USA
Web: <http://www.worldwidediagnostics.com>
E-mail: info@worldwidediagnostics.com

